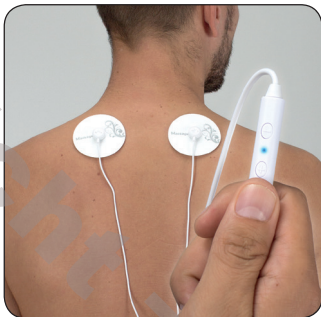
Artikel | Item | Article | Articolo
Artículo | Artikkel | Výrobek | Tuote
Výrobok | Wyrób | продукт | Ürün

4490.5.00

44905GAz_1810

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1 | Introduction

Thank you for deciding for this massager. The technical safety of this product is assured by its compliance to European standards. Read, understand and practice the precautionary and operating instructions.

2 | Designed purpose

It is however not intended to treat medical conditions. This item was conceived for the massage.

- Handy, small, space saving massager
- 6 automatically changing massage modes
- Useable anywhere with smartphone, powerbank or PC as current source
- No charging necessary
- No app necessary

3 | Important information

Retain the manual for future use and forward it to other users!

- Please don't use it in the bathroom or other place with high humidity.
- This stimulator not intended for unattended, personal use by patients who have noncompliant, emotionally disturbed, dementia, or low IQ.
- Children shall not play with the appliance.
- Use of the device does not replace medical consultation and treatment.

3.1 | Precautions

Please do not place pads near heart, genitals or diseased skin.

Do not use in the following cases:

- with use of cardiac defibrillator, pacemaker and diagnosed cardiac problems
- with serious arterial circulatory problems and abnormal blood pressure
- with acute illness thromboses or skin ailments
- with epilepsy, pregnancy, varicose veins, unwellness, acute illnesses and fever
- with numb skin regions
- with chronic pain and reduced pain perception
- do not use electrodes on head
- Stimulation should not be applied over throat or mouth

For additional questions consult your physician!

- Examine the appliance frequently. If there is any damage visible or the supply cord is damaged, it must be replaced by the manufacturer or a qualified person, to avoid a hazard
- The user must not attempt any repairs to the device or any of its accessories
- Manufacturers did not authorize to any maintenance agencies abroad.
- The manufacturer will not be held responsible for the results of maintenance or repairs by unauthorized persons.

4 | Using the device

4.1 | Operating instructions

- a** Connect the two magnetic closures with the gelpads.

 *They can erase information imbedded on magnet strips of EC, credit cards, bank cards, etc.*

- b** Pull the protective film off of the gelpads. stick the gelpads next to each other on the skin area to be massaged
- c** Connect the opposite end of the cable with the power source using the micro-USB plug. Should the plug not fit, the Y-cord is provided as adaptor for USB and USB-C.
- d** The blue LED is on. The device is now ready to be used.

4.2 | Start device

- a** Turn on the device (Button [+])
- b** Press [+] button to adjust the output intensity and the device will go into the work status. There are 8 intensity levels. Press [-] button to decrease the intensity down
- c** Keep press the [-] button to 2 seconds to turn off the device at any mode in the boot state. The blue power LED remains lit until the device is unplugged from the power supply.

4.3 | Automatic shut-off

The appliance automatically switches itself off after 30 minutes.

4.4 | Uses

- There are 6 circulation modes: Massage, Hammer, Relaxation, Scraping, Acupuncture, Cupping
- While the massager is running the 6 massage programs cycle automatically
- The duration of massage should not to exceed 30 minutes. Skin burns may result from using for too long.

- After use, one should wait at least 1 hour before repeating.
- It is recommended not to exceed 1-2 uses per day.

4.5 | Error search

Malfunction	Reason	Countermeasure
Device turns off with use on a power-source	Powersource does not recognize current flow, because the power use of the device is too small therefore shuts off.	Use powersource without automatic use sensing.
Power pulse sensed only weakly or not at all.	Gel pads are used up.	Use new gel pads.

5 | Delivery Contents

Basic device with connection cable	2 gel pads
Y-cable	Instruction manual
Storage box	

6 | Technical data

Model:	C00727	Max. output current:	80mA
Input voltage:	3,3V - 5V	Power:	0,4W
Frequency:	12,8 KHz	Intensity:	1-8 levels
Pulse distance:	200 μ S	Massage modi:	6

7 | Warranty

The guarantee lasts for 24 months from the date of purchase for manufacturers' error in material and quality. If it is used for another purpose or operated incorrectly, no liability for possible damages can be accepted.

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Hydas GmbH & Co.KG
Hydas GmbH & Co. KG
c/o atrikom fulfillment GmbH